

## Invacare® Get-U-Up™ Hydraulic Stand-Up Lift

Comfortable and secure, the Invacare® Get-U-Up™ stand-up lift is the answer to your patient transfer needs. Its ergonomic styling and easy maneuverability make the Get-U-Up an ideal product for every day patient handling. The Get-U-Up helps ensure caregiver security and patient peace of mind. Ideal for use with partial weight-bearing patients as well as those needing rehabilitation support, the Invacare Get-U-Up stand-up lift offers safety, comfort and stability.



Yes, you can.®




Model No. GHS350

### Features


- Innovative adjustment features adapt to a wide variety of body sizes and shapes. That means comfort at any height, width and weight range up to 350 lb.
- Low friction casters significantly improve rollability and contribute to caregiver safety
- Multi-function slings allow for quick toileting, stand-assist or full-support seated transfer
- Adjustable leg support features a soft pad with lateral contour that secures legs in multiple positions for comfort at any height.
- Retractable, non-slip footplate provides stable base for transfer
- Under-bed clearance of 4.5" allows accessibility to most homecare style beds
- Ergonomic features make it ideal for toileting transfers, as well as dependent and rehabilitation purposes
- Fast, simple sling attachment that protects against accidental disengagement; wide push handle, good access for clothing removal; and comfortable slings with extended back support that prevents pressure in the under arm area
- Pump handle can rotate from side to side for convenience of caregiver
- Limited warranty includes three years on the lift, one year on the hydraulic pump

### Specifications

Height	
Maximum at sling hookup	65"
Minimum at sling hookup	36"
Base height (clearance)	4.5"
Base width	
Open	37"
Closed	23"
Base length	
	40.5"
Casters	
Front	3"
Rear	5"
Weight capacity	
	350 lb.
Product weight	
Out of carton	88 lb.
Shipping	100 lb.

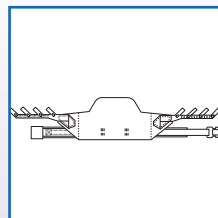

Model no. R130 Standing sling.

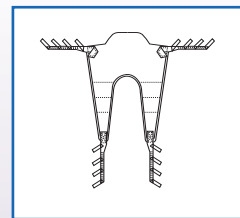

Model no. R134 Transfer sling.

### Sling Specifications

	R130	R134
Width	44"	34"
Length	16"	42"
Product weight		
Out of carton	1.1 lb.	1.3 lb.
Shipping	1.6 lb.	1.8 lb.

Invacare Corporation  
www.invacare.com

### USA

One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900


### Canada

570 Matheson Blvd. E.  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

All rights reserved.  
Trademarks are identified by the symbols ™ and ®.  
All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted.

© 2007  
Invacare Corporation

Form No. 07-075

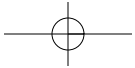

Invacare® Get-U-Up™ Hydraulic Stand-Up Lift

Personal Care • Manual Wheelchairs  
HOMECARE  
Bed • TSS • Ambulifters •

# Get-U-Up™



Model No. GHS350

