



INVACARE® MICROAIR® MA1000
THERAPEUTIC SUPPORT SURFACE



Yes, you can.®

microAIR®

Invacare® microAIR® MA1000

Alternating Pressure Low Air Loss Mattress System

Medium - High Risk
Stage I-IV

3-in-1 Alternating Pressure System combines True Low Air Loss and advanced clinical features, creating the Ultimate Solution for Patient Comfort and Pressure Relief



Highly vapor permeable and oversized pliable quilted nylon top cover provides low shear, friction and moisture protection

Added convenience of **Integrated side bolsters** that can be inflated and deflated on-demand




Nylon non-slip base cover helps reduce shear and friction during bed articulation



Advanced 3-in-1 alternating function inflates and deflates cells in a 3:1 cycle, providing active pressure relief by supporting 2/3 of the body at any given time

Features & Benefits

10 adjustable comfort settings depending on the patient's therapy requirement

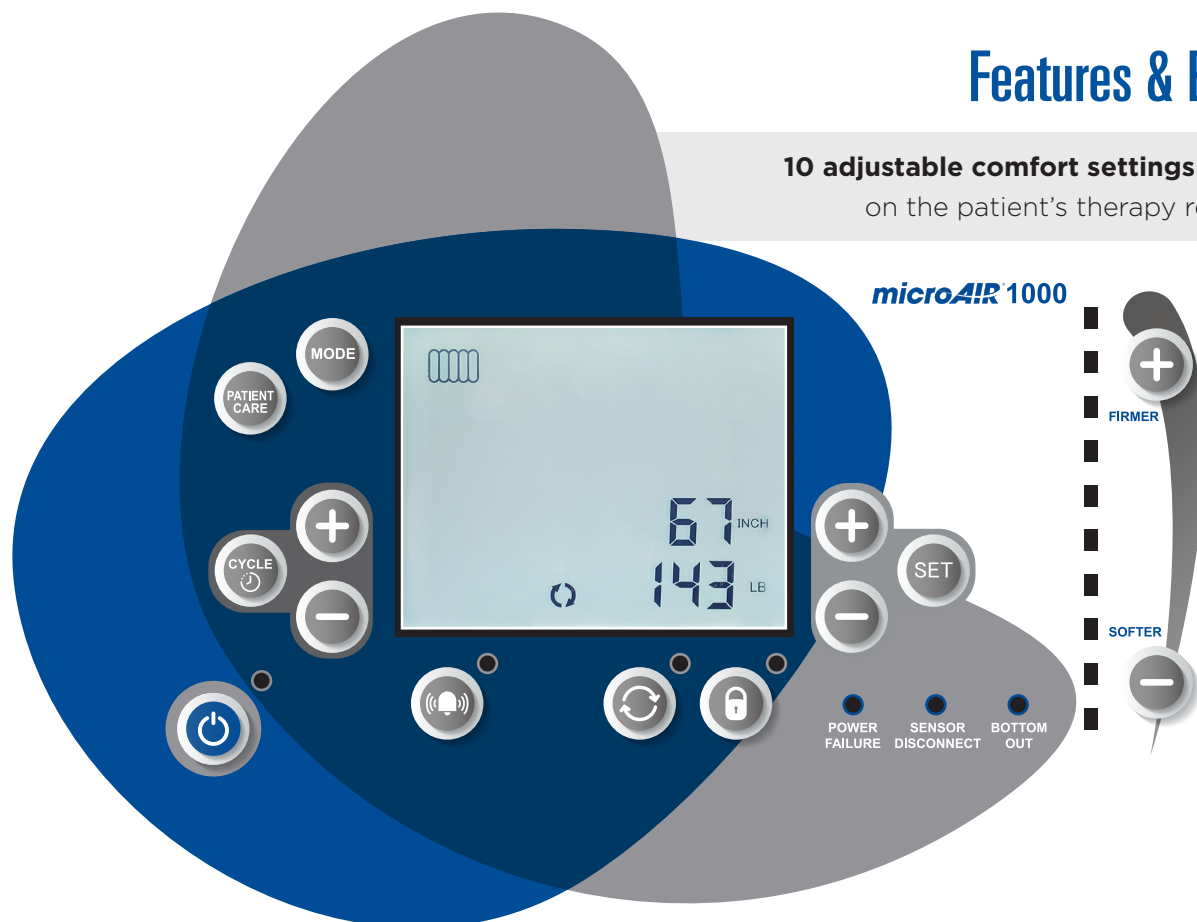

Auto Set Mode automatically sets optimal pressure based on patient's weight, height and position for constant pressure monitoring

60 second pulsation feature helps improve capillary blood flow by continually readjusting to a level surface pressure.

Auto Fowler Mode safety feature can detect the pressure of patient's buttocks automatically and adjust pressure accordingly

True Low Air Loss Therapy provides extra comfort for moisture and temperature control

Patient Care Mode provides quick maximum inflation within seconds to aid with transfers and nursing procedures

Customizable cycle times to meet patient's varying needs



Foam base helps prevent bottoming out in the event of a power failure

Soft nylon urethane backed cells provide easy cleaning and helps minimize the proliferation of bacterial growth

HeelSense™ Technology helps provide further therapy and comfort by decreasing pressure in patient's vulnerable heel area

Request a Quote or Place an Order Today!



Contact your Inside Sales Representative

or




Contact Customer Service at 800.333.6900

Invacare Corporation
www.invacare.com

USA
One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada
570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

Invacare microAIR MA1000 Alternating Pressure Low Air Loss Mattress System

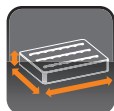

Model #	Dimensions	Weight Capacity	Flow Rate (Pump)	Therapy Type	Alarms	Limited Warranty
MA1000	80" L 36" W 10" H	600 lb.	1275 LPM	Alternating Pressure Pulsation True Low Air Loss	Audio & Visual	1 Year (Cover/Mattress) 1 Years (Power Unit)
MA1000M	microAIR MA1000 Mattress only					
MA1000P	microAIR MA1000 Pump only					

©2019
Invacare Corporation.
All rights reserved.
Trademarks are identified by the symbols™ and ©. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.
Form No. 19-030
190209



Yes, you can.®