

INVACARE® WALKERS



With the **Invacare® Walker**, providers, clinicians and consumers gain the stability they need at an economical price. Invacare Walkers offer a wide, deep frame with a large number of height adjustments and are offered in junior and adult sizes. The walker has a lower side brace for added stability and is stable, lightweight and easy to lift and maneuver. Walkers include PVC handgrips that provide comfort and long-lasting wear and anti-rattle silencers for quiet operation.

Invacare® Dual Blue-Release™ Adult Walker Features

Lightweight, wide, deep aluminum frame

Dual Blue-Release mechanisms provide both visual and audible “locked” cues

Rigid dual side brace design

Large number of height adjustments

Anti-rattle “silencers” for quiet operation

PVC handgrips for comfort and long lasting wear

Compatible with most Invacare Walker Accessories (not to be used with Walker Tray or Walker Pouch)

HCPCS Code:

E0135 Adult/Jr.

E0143 Wheeled Adult/Jr.



1. Dual Blue-Release button



6240-A: Dual Blue-Release Adult Walker



6240-5F: Dual Blue-Release Adult Walker with 5" Wheels

Invacare® Single-Release Walker Features

Lightweight aluminum frame

One button to press and fold

Single button provides both visual and audible “locked” cues

Large range of height adjustments

PVC handgrips for comfort and long-lasting wear

Compatible with most Invacare Walker Accessories (not to be used with Walker Tray or Walker Pouch)

HCPCS Code:

E0135 Adult/Jr.

E0143 Wheeled Adult/Jr.



1. Single-Release button



6281-A: Single-Release Paddle Adult Walker



6281-JR: Single-Release Paddle Jr. Walker

Invacare® I-Class™ Paddle Walker Features

Composite lower side brace stiffens walker frame and adds stability

Folding mechanisms are easy to reach from a seated position

Lighter and easier to lift and maneuver


PVS handgrips wear longer and provide greater comfort

Walkers fold easily for transport or storage

HCPCS Code:

E0135 Adult/Jr./Single

E0143 Wheeled Adult/Jr.

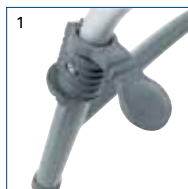

1. Dual-Release Paddle button



6291-A: Dual-Release Paddle Adult Walker



6291-5F: Dual-Release Adult Paddle Walker with 5" Wheels

The **Invacare® Heavy-Duty Dual-Release Paddle Walker** is strong, stable and simple to use, helping individuals to walk with confidence and ease. This walker also is compatible with Invacare leg accessories.

Features


Compatible with model # 6375 bariatric 5" fixed wheels and leg extensions

The **Invacare® Bariatric Adult Walker** provides a wide and deep frame that can support an individual up to 700 lb. The dual-release buttons provide both visual and audible "locked" cues for the user. The durable aluminum frame provides maximum stability for consumers and the legs offer variable height adjustments.

Features

Compatible with Invacare model # 6372 bariatric 5" fixed wheels and leg extensions

HCPCS Code: E0148

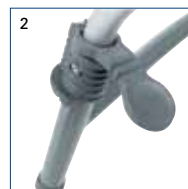

1. Dual Blue-Release button (model 6291-HDA)

2. Dual-Release Paddle button (model 6441-A)

3. Bariatric size



6291-HDA: Heavy-Duty Dual-Release Adult Walker



6441-A: Bariatric Dual Release Adult Walker

The **Invacare® Two Step Walker** features a wide, deep frame with a large number of height adjustments and a lower side brace for added stability. It features two stages of grips to help the user ease from a sitting to standing position. The Invacare Two Step Walker is light, easy to lift and maneuver, and features foam hand grips for comfort and anti-rattle silencers for quiet operation.

Features

Combination single-release folding walker and sit to stand assistance device

Can also be used as a toilet safety frame

Adjustable height settings

Accessories currently not available

HCPCS Code: E0135



The **Invacare® Adult Hemi Walker** is designed for patients with limited or no dexterity in one arm or hand. This lightweight walker combines support and stability, yet is easy to lift and maneuver, even with one hand.

Features

Wider base provides more stability than a quad cane

Anodized aluminum tubing, lightweight and maneuverable

Compact size for easy storage and transport

Extra handgrip provides ease in folding

Slip-resistant rubber tips

HCPCS Code: E0135



Walker Specifications



Model #	Width Inside Hand Grips	Width Inside Base Legs	Depth Opened At Base	Depth Folded	Height Adjustment	Product Weight	Weight Capacity
6240-A	17"	21"	16"	4"	32"-39"	5 lb. 7 oz.	300 lb.
6240-JR	17"	20"	16"	4"	25"-32"	5 lb. 5 oz.	300 lb.
6240-5F	17"	21"	18.5"	4"	32"-39"	6 lb. 6 oz.	300 lb.
6240-JR5F	17"	20"	18"	4"	25"-32"	6 lb. 12 oz.	300 lb.
6281-A	17"	21"	17.5"	4.5"	30.375"-37.375"	6 lb. 2 oz.	300 lb.
6281-JR	17"	21.5"	16.5"	4.5"	25.375"-32.375"	5 lb. 14 oz.	300 lb.
6291-A	17"	21.5"	17.5"	4"	30.375"-37.375"	5 lb. 5 oz.	300 lb.
6291-1	17"	21.5"	17.5"	4"	30.375"-37.375"	5 lb. 5 oz.	300 lb.
6291-JR	17"	20.5"	16.5"	4"	25.375"-32.375"	5 lb. 1 oz.	300 lb.
6291-3F	17"	21.5"	17.5"	4.5"	32"-38"	5 lb. 5oz.	300 lb.
6291-JR3F	17"	20.25"	16.5"	4.5"	27"-33"	5 lb. 1oz.	300 lb.
6291-5F**	17"	21.5"	17.5"	4.5"	33"-39"	5 lb. 6 oz.	300 lb.
6291-JR5F**	17"	20.25"	16.5"	4.5"	28"-34"	5 lb. 6 oz.	300 lb.
6291-HDA	20"	24"	17.5"	4"	30.375"-37.375"	6 lb. 11 oz.	500 lb.
6441-A	20"	23.5"	16"	4.25"	32.25"-41.25"	7 lb.	700 lb.
P421B	16.5"	21.5"	19"	5"	31"-36"	5 lb. 8 oz.	286 lb.
6252	n/a	n/a	n/a	3"	29.5"-34.5"	2 lb. 14 oz.	250 lb.



	Patient Height Range*	Units Per Carton	Limited Warranty
6240-A	5'4"-6'6"	4	5 Years
6240-JR	4'4"-5'7"	4	5 Years
6240-5F	5'4"-6'6"	4	5 Years
6240-JR5F	4'4"-5'7"	4	5 Years
6281-A	5'3"-6'4"	4	5 Years
6281-JR	4'4"-5'7"	4	5 Years
6291-A	5'3"-6'4"	4	5 Years
6291-1	5'3"-6'4"	1	5 Years
6291-JR	4'4"-5'7"	4	5 Years
6291-3F	5'4"-6'4"	4	5 Years
6291-JR3F	4'6"-5'6"	4	5 Years
6291-5F**	5'6"-6'6"	4	5 Years
6291-JR5F**	4'4"-5'6"	4	5 Years
6291-HDA	5'3"-6'4"	2	5 Years
6441-A	5'3"-6'10"	2	5 Years
P421B	5'3"-6'4"	5	5 Years
6252	5'1"-5'10"	2	5 Years

* All patient height ranges are approximate. Each individual should always consult with his/her physician or therapist to determine proper adjustment and usage.

** With 13" rear leg extensions.

"F" designation in Model no. 3" or 5" fixed wheels with glide tips.



1. **Walker Basket**

Model # 6096

For model # 6240-A, 6240-JR, 6240-5F, 6240-JR5F walkers only

2. **Walker Basket**

Model # 6098

Compatible with all walkers except model # 6291-HDA, 6252, 6300-JRA, 6300-ATA

3. **Walker Tray**

Model # 6007

Compatible with all walkers except model # 6291-HDA, 6252, 6300-JRA, 6300-ATA, 6240, 6281 walkers

4. **Walker Tray**

Model # 6002

For model # 6240-A, 6240-JR, 6240-5F, 6240-JR5F walkers only

5. **Walker Pouch**

Model # 6015-4

Compatible with all walkers except model # 6291-HDA, 6252, 6300-JRA, 6300-ATA

6. **Platform Attachment**

Model # 6027

Compatible with all walkers except model # 6291-HDA, 6441-A, 6252 or 6300-JRA, 6300-ATA



1. **Invacare® Court-Side Glide™ Attachment**

Model # 6096

For model # 6240-A, 6240-JR, 6240-5F, 6240-JR5F walkers only




1. **Glide Tips**

Model # 6280

Packaged 24 per carton

Durable construction for long wear

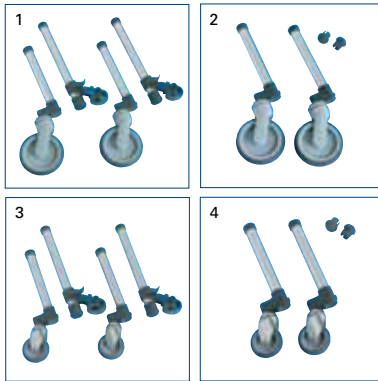


1. **5" Swivel Wheels with 13" Cambered Invacare® SureGlide™ Brakes**
Model # 6268
2. **5" Swivel Wheel Attachments with Rear Glide Tips**
Model # 6267
3. **3" Swivel Wheels with 13" Cambered Invacare® SureGlide™ Brakes**
Model # 6272
4. **3" Swivel Wheel Attachments with Rear Glide Tips**
Model # 6273




**Also available in retail packaging with a minimum purchase of four sets. Order as Model # 6270-4 and 6271-4.

1. **3" Single Fixed-Wheel Attachments with Rear Glide Tips**
Model # 6270**
2. **5" Single Fixed Wheel Attachments with Rear Glide Tips**
Model # 6271**
Use 13" rear leg extensions for maximum adjustability
3. **Tall 3" Single Fixed-Wheel Attachments with Tall Rear Legs and Glide Tips**
Model # 6278 with 13" rear leg extensions
Adjusts to accommodate persons up to 6'6" tall
4. **Bariatric Wheel Kit**
Model # 6372 for use with 6441-A. 6375 for 6291-HDA not shown
Dual 5" fixed wheels for added stability
Large number of height adjustments

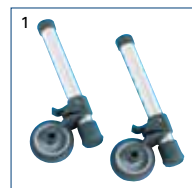

1. **10" Rear Wheel Brakes**
Model # 6264
Use with 3" or 5" single fixed-wheel attachments
2. **13" Rear Wheel Brakes**
Model # 6265 (not pictured)
Use with 3" or 5" swivel wheel attachments or tall 3" single fixed-wheel attachments

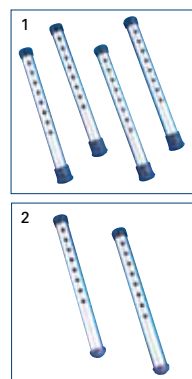

1. **Tall Leg Extensions**
Model # 6275
Adjusts to accommodate persons up to 6' 6" tall
2. **13" Rear Leg Extensions with Glide Tips**
Model # 6266
Use with swivel or fixed wheels

