

Compact RWD Powerbase Wheelchair

P310 Regal

The unique compact size is designed for maximum maneuverability in narrow spaces and walkways. Ideal powerchair for indoor and outdoor use.

Technical Specs

max.speed : up to 5mph
 range : up to 18mi.
 turning radius : 21"
 ground clearance : 3"
 motor : DC 24V, 180W
 controller : PG VR2 50A
 battery : 12V / U1(34AH)x2
 charger : 5A off-board
 gradient : 12 °
 caster : 8" foam filled tires
 drive wheel : 10" foam filled tires
 brake : intelligent, regenerative, electromagnetic brakes

Features

- Stylish and compact rear wheel drive model with 10" drive wheels and alloy rims.
- Extremely low center of gravity.
- Lower seat height which facilitates transfers and makes it the most stable unit in its class.
- Sturdy, deluxe captain seat with angle adjustable back and flip-up armrests.
- Excellent maneuverability.
- Programmable controller with a wide range of selections to suit the needs and demands of different users.
- Off-board charger.
- Semi-reclining seat.
- Colors available : Red, Blue, Black.

Options

Elevating leg rests. Oxygen tank holder.
 Cane & crutch holder. Walker holder.
 Rear foldable basket. Canopy. Cup holder. Side bag.




Weight (lb)



batteries	capacity
w/o 137	300
w/ 189	

Dimensions (in)

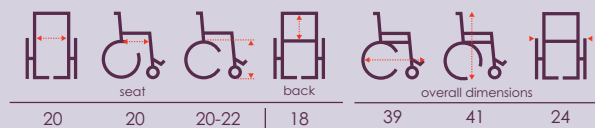

Medicare Code

K0823 P310-2S C
 K0822 P310-2S S