

# Protekt™ Aire 9000

Low Air Loss and Pulsating  
Mattress System




## FEATURES AND BENEFITS:

- 1275 LPM (liters per minute) true low air loss.
- Provides effective true low air loss relief, or pulse therapy along with continuous low air loss relief.
- In pulse mode every minute the mattress pressure is dropped to 80% of the set pressure for few seconds.
- Max Flow mode with 15 minute timer for rapid inflation of the mattress, inflates in 15 ~ 50 seconds depending on the size of the mattress.
- In upright (fowler) mode the pump inflates the mattress to 50% higher than the set pressure, to eliminate patient bottoming.
- Lock out key completely locks the control panel to eliminate tampering of the patient settings, including the power switch.
- Audio and visual power failure alarms.
- User-friendly control panel with dome keys and 7 segment LED display for simple user interface, displays current mode, and patient comfort pressure level from Soft 1 to Firm 10 ( $8 \pm 4$  to  $32 \pm 6$  mmHg).
- Low friction and low shear force producing, urethane coated breathable nylon top sheet, which is liquid proof and highly vapor permeable.
- 2" convoluted safety foam base to provide additional support to the patient in the event of power failure.
- Weight capacity: 1,000 lbs.
- Mattress dimensions: 36"x80"x10".
- UL/CSA approved power cord.
- 18 month warranty.

80090	Protekt™ Aire 9000 BL
-------	-----------------------

# Protekt O2 Cushion

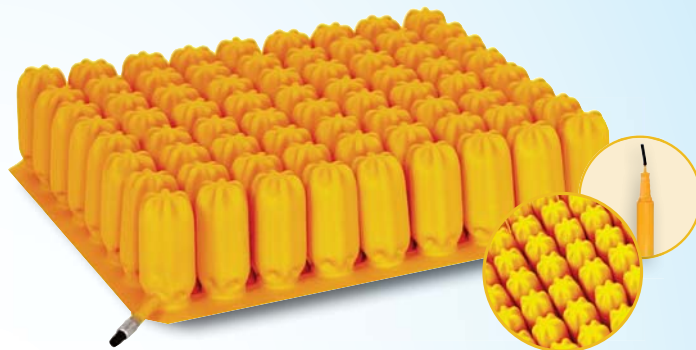

## FEATURES AND BENEFITS:

The Protekt O2 cushion is Proactive Medicals premier wheelchair cushion for the prevention and treatment of pressure sores. The soft interconnected air cells allow for deep immersion providing constant low pressure across the seating surface.

- Air adjustable 4" or 2" cells offer deep immersion to provide excellent pressure redistribution for tissue protection and treatment.
- Dual valve compartments address seating challenges such as pelvic obliquities, sacral seating and enhances posture control.
- Dual compartments can be adjusted independently to increase positioning and stability for side to side control.
- Nylon breathable cover with non-skid bottom and safety straps.
- Easy to use hand pump.
- Flame resistant, latex free neoprene.
- 2 year warranty.
- Patch kit included.

78001	16"x16"x2"
78002	18"x16"x2"
78003	20"x16"x2"
78004	16"x16"x4"
78005	18"x16"x4"
78006	20"x16"x4"