


Geo-Mattress® UltraMax™



Now made with

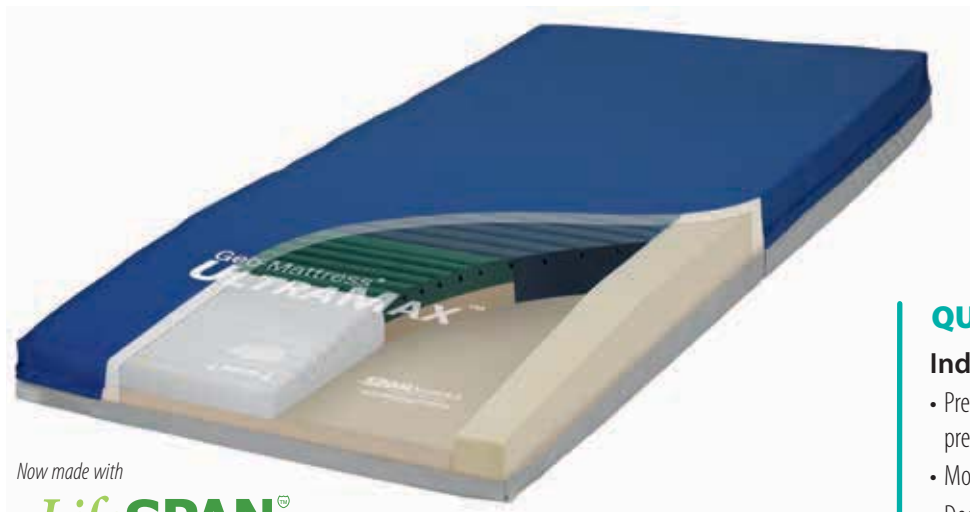

bleach-resistant fabric cover
(see page 15 for details)



Designed to exceed even the most rigorous requirements of the acute care environment.

Precision-engineered to deliver critical performance advantages inside and out! On the outside, the individually articulating cells of the clinically proven Geo-Matt® top surface combine with Span-America's proprietary Shear Transfer Zones® design to provide unmatched protection from the harmful effects of micro shear, macro shear and rotational shear. Helps reduce caregiver injury by "gliding" the user back to position following head of bed elevation. On the inside, the 3-Dimensional Zoning™ core design of the UltraMax provides anatomically calibrated levels of support and unmatched user comfort. And in the seat area – where the mechanical fatigue of bed gatching and concentrated user weight affect durability most – the unique cell structure of the high resiliency UHP 500™ layer provides a 55% higher Support Factor than found in typical medical grade foams. Viscoelastic foam and gentle Heel Slope® provide added protection at vulnerable heels.

Therapy Provided

- ✓ Non-powered Pressure Redistribution

QUICK FACTS

Indications:

- Prevention and early intervention of Stage 1-4 pressure injuries¹
- Moderate to high risk for skin breakdown
- Decreased mobility
- Repositions frequently by self or caregiver

Key Advantages:

- Minimized shearing thru Geo-Matt® surface geometry and Shear Transfer Zones® within stretch cover (see description on page 15)
- Outstanding comfort for longer, more restful sleep
- Superior in-bed mobility, caregiver ergonomics
- Enhanced durability, protection from "bottoming out"

Item #	Description
UMX7536-29	75"L x 36"W x 6"H
UMX8036-29	80"L x 36"W x 6"H
UMX8436-29	84"L x 36"W x 6"H
UMX85V-29	For recessed deck frames (see page 19)
UMX80CA-29	80" L for CareAssist®+
UMX84CA-29	84" L for CareAssist®+
* UMX8039-29	80"L x 39"W x 6"H
* UMX8042-29	80"L x 42"W x 6"H
* UMX8439-29	84"L x 39"W x 6"H
* UMX8442-29	84"L x 42"W x 6"H
* Allow up to 10 working days lead time for these sizes.	
Weight:	24 lbs.
Mattress Weight Limit:	500 lbs (36"); 600 lbs (39" and 42")
Warranty:	8 years, all components. Non-prorated
Flammability:	All models conform to NFPA 101 (Life Safety Code™); CPSC 16 CFR Parts 1632 and 1633; Cal TB 117 & 129

¹ May be appropriate for multiple stage 2 or single stage 3 injuries, based on full assessment of skin status and repositioning required, according to best clinical practice and judgment

+ Hill-Rom is a mark, and CareAssist® is a registered mark, of Hill-Rom® Services, Inc.





Geo-Mattress® UltraMax™ for Stretchers

Category:
Premium Pressure Redistribution
Support Surface




3-Dimensional Zoning™ design provides anatomically-calibrated levels of support and immersion throughout.



Extra-soft heel section facilitates immersion to provide added heel protection.

Design Advantages Inside...

Engineered, multi-layer
support core

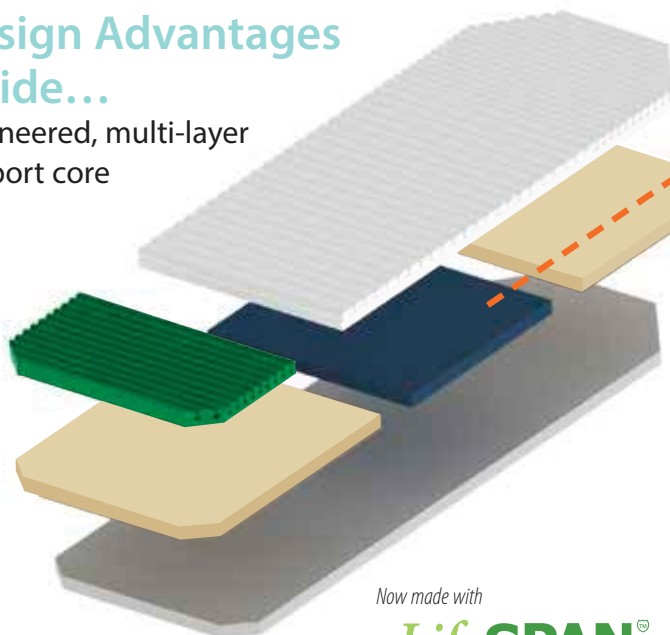

UHP 500™ Ultra High Performance foam in high-fatigue seat section; 55% higher support factor yields superior resiliency and recovery.



Geo-Matt® non-shearing surface creates individually articulating cells; Ring-of-Air® channels dissipate heat and moisture.

Now made with



bleach-resistant fabric cover
(see page 15 for details)

Specifications

[All models listed are 5" high to facilitate maximum immersion]

Item #	Dimensions, overall (L)ength and (W)idth	(H)ead-end taper length	(F)oot-end taper length
47318	75"L x 26"W	7"	0"
49297	75"L x 26"W	13"	4"
97107	76"L x 26"W	10"	4"
10635	76"L x 26"W	12"	5"
2879529	76"L x 27"W	11"	5"
46234	76"L x 29"W	10"	8"
79012	76"L x 30"W	13"	4"
98969	78"L x 27"W	12"	12"
81921	2-pc "Gynie"	n/a	n/a

Other dimensions also available upon request.


Supports: up to 700 lbs.

Warranty: 3-year, non-prorated

Mattress Weight: approx. 20 lbs.

All components: 100% latex-free and 100% PBDE-free

Flammability: Conforms to NFPA 101 (Life Safety Code™); CPSC 16 CFR Parts 1632 and 1633; Cal TB 117 & 129.

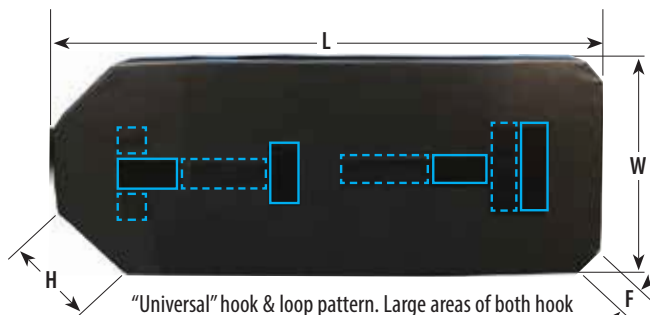

"Universal" hook & loop pattern. Large areas of both hook and loop on bottom of mattress help in anchoring it to a wide variety of stretcher models.

