

PressureGuard® APM²

PressureGuard® APM² with digital control unit



Now made with
LifeSPAN™
fortified healthcare fabric
bleach-resistant fabric cover
(see page 15 for details)

¹ In situations where positioning on the pressure ulcer cannot be avoided (e.g. the individual has multiple ulcers on multiple surfaces), limit the amount of time the individual is positioned on the ulcer. Source: NPAUP/EPLUAP Pressure Ulcer Prevention and Treatment, Clinical Practice Guideline, p 66-67

Therapy Provided

- ✓ Powered Pressure Redistribution
- ✓ Alternating Pressure
- ✓ Lateral Rotation

The APM² has been an industry-leading air therapy solution for more than a decade, based in large part on its documented superiority vs. typical, roll-up treatment surfaces. Newly upgraded cover design and control unit provide even greater immersion, envelopment, shear reduction, treatment flexibility, caregiver ergonomics and user comfort.


QUICK FACTS

Indications:

- Treatment of Stage 1-4 pressure injuries
- Patients that can't/won't reposition frequently
- Can be positioned off pressure injury in at least two positions¹
- Prevention for high risk patients

– New, Upgraded Features –

Digital Multi-function Control Unit

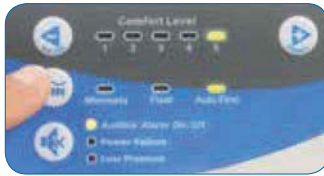

Lightweight and quiet, it delivers an unmatched combination of treatment flexibility, ease of use, and dependability.

- Low Pressure indicator light and selectable audible alarm.
- Power Failure indicator light and selectable audible alarm.
- Simple, indexed comfort adjustment.
- Timed (20-minute) "Auto Firm" mode for transfers, ADLs, & CPR.
- "Float" mode provides powered flotation therapy.
- "Alternate" mode changes loading in 10-minute cycle; controlled at mattress for choice of lateral rotation or alternating pressure therapy.




Shear Transfer Zones® Cover Design (US Patent Number 8,438,682)

In addition to the Geo-Matt® segmented top surface, Span-America's exclusive Shear Transfer Zones® provides an additional measure of shearing protection in the form of the silicone-coated, shear-minimizing fabric bands located on the underside of the bi-directional stretch cover.

This Shear Transfer Zones design helps:

- Prevent heels, sacrum, and scapula from "digging into" surface.
- "Glide" user back to original position following HOB elevation.
- Protect against damaging effects of micro shear, macro shear, and rotational (pivot-induced) shear.
- Provide patient stability by transferring shear to more shear-tolerant anchor points.

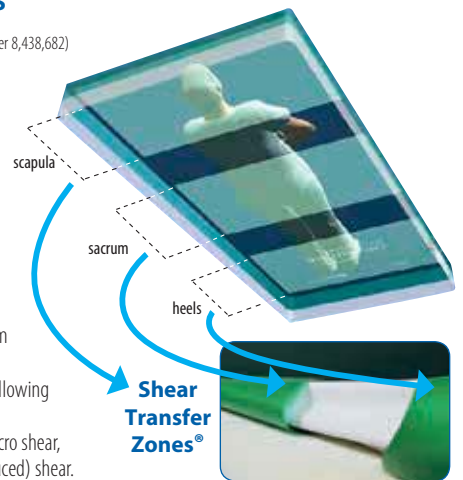

Item #	Description	Dimensions
5875LR-29	APM ² with digital control unit	75"L x 35"W x 7"H
5880LR-29	APM ² with digital control unit	80"L x 35"W x 7"H
5884LR-29	APM ² with digital control unit	84"L x 35"W x 7"H
5900	Digital control unit only	
Control Unit: Weight: 6.3 lbs.		Dimensions: 12" x 7.5" x 5"
Voltage: 120 AC		Max current: 1.0 Amp
Frequency: 60 Hertz		UL Listing: Certified to IEC 60601-1-1

Medicare Approved: E0277 (Group II)
Weight limit: 350 lbs
Mattress weight: approx. 20 lbs.
Warranty: 18 months (all components)

All models are fire barrier version (FBV) as indicated by "-29" suffix, denoting certified conformance to NFPA 101 (Life Safety Code); Cal Tech Bulletin #117; 16 CFR Parts 1632 and 1633; and BFD IX-11.



Extra-wide and Consumer Sizes

Several Span models are available in sizes wider than typical healthcare beds. These non-bariatric options provide comfortable, effective use by those who require a wider surface but do not require a bariatric weight capacity. Options include:

PressureGuard® powered air therapy (treatment):

Easy Air®

Supports up to 500 lbs.
Widths: 39"
See page 4



Easy Air® LR

Supports up to 500 lbs.
Widths: 39", 53"
See page 5



Protocol®

Supports up to 500 lbs.
Widths: 42"
See page 6



Topper®

Supports up to 1000 lbs.
Widths: 42", 48", 54", 60"
See page 7



PressureGuard® convertible and non-powered air therapy (treatment):

Custom Care® Convertible

Supports up to 500 lbs.
Widths: 42"
See page 16



Span Care® Convertible

Supports up to 700 lbs.
Widths: 42"
See page 13



Custom Care®

Supports up to 500 lbs.
Widths: 42"
See page 18



CFT®


For users 100 – 750 lbs.
Widths: 42", 48", 54"
See page 12



Geo-Mattress® therapeutic foam (prevention):


UltraMax™

Supports up to 600 lbs.
Widths: 39", 42"
See page 26

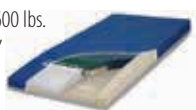

Atlas™

Supports up to 750 lbs.
Widths: 39", 42", 48", 54"
See page 29

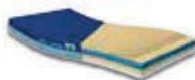

Max®, Plus® and Pro™ models

Supports up to 750 lbs.
Widths: 39", 42", 48", 54"
See page 28



Also available....PressureGuard In-Home Styles. Extra-wide, consumer-sized models accommodate one user or two.



Therapy on Patient Right

When ordering, customer must specify if the dynamic air therapy should be positioned in the **Center** of the mattress (for single user), or on the **Patient's Right Side** (for two users or a single user who cannot/prefers not to move to the center of the mattress). Note: "Patient Right" is with patient in a backlying position.



Therapy in Center

PressureGuard® APM² with digital control unit (see page 8 for product description)

Size/Dimensions	Therapy on Patient Right	Therapy in Center	Weight limit
Dual King* / 80"L x 38"W	AP28038-29		350 lbs. each user
Full / 75"L x 54"W	AP27554-RT-29	AP27554-CTR-29	350 lbs. each user
Full Long / 80"L x 54"W	AP28054-RT-29	AP28054-CTR-29	350 lbs. each user
Queen / 80"L x 60"W	AP28060-RT-29	AP28060-CTR-29	350 lbs. each user

PressureGuard CFT® (see page 12 for product description)

Size/Dimensions	Therapy across entire surface	Weight limit
Dual King* / 80"L x 38"W	CF8038-HS-29	500 lbs.
Full / 75"L x 54"W	CF7554-HS-29	750 lbs. total**
Full Long / 80"L x 54"W	CF8054-HS-29	750 lbs. total**
Queen / 80"L x 60"W	CF8060-HS-29	750 lbs. total**

* To accommodate a Dual (split) King surface, 80"L x 76" W, order this item together with our matching non-powered air mattress (RW8038-29) for the other half of the bed.

** 1 or 2 users.

Non-prorated warranty: 18 months (APM² Deluxe); 5 years (CFT).

Cover cleaning: Wipe clean with neutral suds or approved cleaning agents. Avoid bleach.

Lead times: Allow up to 10 working days lead time for all sizes.

Shipping: All In-Home Style units ship fully assembled in a flat carton via common carrier.

UL Listing: APM² Deluxe conforms to UL 544.

See individual product brochures for complete technical specifications.

Need
Microclimate
Management
("Low Air Loss")?
Add the Topper!
(See page 7)

